HappyLight

Model: VT20

Edit#: 002	Title: VT20 Manual Trim to A5
Date: 09/06/17	Version: Rev20

Dear Customer,

Thank you for purchasing the Verilux® HappyLight® Liberty™ Series Energy Lamp. You now own an innovative product, manufactured to the highest standards and backed by a one-year limited warranty.

This product brings Natural Spectrum Daylight into your home or office. Many other healthy-lighting products are available online. Visit us on the web at www.verilux.com to learn more about all of our quality Verilux products, or call our toll-free number: 1-800-786-6850.

As a Verilux customer, your satisfaction means everything to us. We look forward to serving you now and in the future.

Have a Bright Day!

Nicholas Harmon
President
Verilux, Inc.

Important Safeguards 4
Features. 5
Benefits 6
Device Pretest 6
Operation. 7
 Selecting a Location

Edit#: 002	Title: VT20 Manual Trim to A5
Date: 09/06/17	Version: Rev20

- To avoid electrocution, do not operate this device near water.

▲WARNING:

- Do not use with power supply voltage other than 120 VAC.
- Do not defeat the safety purpose of the polarized plug. A polarized plug has two blades with one wider than the other. The wider blade is provided for your safety. If the provided plug does not fit in your outlet, consult an electrician for replacement of the obsolete outlet.
- Close supervision is recommended when these lamps are used by or near children, or others who may not understand the need for these precautions.
- Do not operate lamp in close proximity to flammable or combustible vapors, such as aerosol spray products, or where oxygen is being administered.
- To prevent risk of shock or personal injury when cleaning lamp, make sure you have turned off and unplugged the lamp and allowed sufficient time for the bulb to cool.
- Do not block vents. Vents are necessary for proper function of the electronic ballast; blocking them may cause overheating of the electronics and increase risk of fire.
- Do not shorten or cut the power cord.
- Always unplug lamp before performing bulb replacement or installation and allow sufficient time for the bulb to cool.
- Light bulbs get hot quickly! Do not touch the bulb while the lamp is on.
- This bulb contains a trace amount of mercury. For safe handling, cleanup, disposal and recycling in accordance with applicable laws in US see 111.1, a. /c

▲ CAUTION: — continued

- This lamp may cause interference with radios, cordless telephones or devices that use a wireless remote control, such as televisions. If interference occurs, move the lamp away from the device, plug the lamp or the device into a different outlet or move the lamp out of the line of sight of the remote control receiver.
- Unplug lamp during lightning storms or when unused for long periods.
-

Features

Edit#: 002	Title: VT20 Manual Trim to A5
Date: 09/06/17	Version: Rev20

Verilux HappyLights provide

Device Pretest

To help determine if you have sensitivity to light treatment, we recommend that you perform this pretest upon first use of this device:

- Step 1 - Start this test with the comfort lens. Turn on the lamp and position it about 18 inches from your face. (Do not look directly into the lamp.)
- Step 2 - Operate the lamp for approximately 10 minutes and turn it off. If the light has caused any discomfort, stop using the lamp and consult a doctor. If you have no discomfort, you may consider gradually decreasing the distance from the lamp and increasing the use time, which results in more lamp intensity and greater effect. You also may change from the comfort lens to the high-energy lens and repeat the device pretest. For your initial sessions, you should monitor the degree of any discomfort and accordingly adjust the time used, lens used and lamp location.

▲WARNING: Some people may experience light sensitivity to fluorescent lighting. Discontinue use if you experience any discomfort.

▲WARNING: Certain medications may induce light sensitivity (e.g., anti-depressants) which may result in discomfort or a level of hyperactivity. Consult with your health care provider before using this product.

▲WARNING: If you have a pre-existing eye condition, contact your health care provider prior to use.

Selecting a Location

- Place your Verilux HappyLight Energy Lamp on a table or other flat surface above, at or below eye level depending on your desired result.
-

Instructions for Use

Use consists of sitting close to the HappyLight Energy Lamp with your eyes open. While staying oriented toward the light, feel free to engage in other activities such as reading, writing, hobbies and crafts, sewing, needlepoint, applying make-up or working at your computer or desk. You may tilt the lamp as required.

Edit#: 002	Title: VT20 Manual Trim to A5
Date: 09/06/17	Version: Rev20

Your HappyLight uses a bulb manufactured only by Verilux. To order a replacement bulb, please call Verilux, toll-free, at 1-800-454-4408 or order online at <http://www.verilux.com/customer-service>

To Replace the Bulb:

1.


4.

Edit#: 002	Title: VT20 Manual Trim to A5
Date: 09/06/17	Version: Rev20

Changing the Lenses/Lens Replacement

Your HappyLight Liberty Series Energy Lamp comes with interchangeable comfort and high-energy lenses, inserted in a bezel and held with magnets to the lamp body, that regulate the light intensity, reduce glare and provide your eyes with even stimulation. You can choose – or change – the covers to get the amount of light you want for your light therapy program:

1. Remove the lens by inserting a coin or your finger into the opening between the lens bezel and the lamp body (A).
2. Gently pull the lens and bezel away from the lamp body, lifting slowly as you pull (B).
3. Hold the bezel at the top and pull away the lens (C).
4. To install the other lens, work from the underside of the bezel and lens. Align the six cutouts in the lens (D) with the six magnet clips (E). Lower the lens, and slide it back into position.

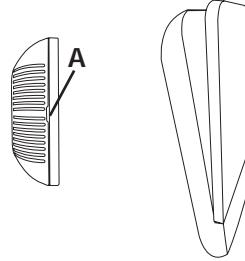⚠ WARNING: This product contains magnets. If you have a pacemaker, implantable defibrillator, insulin pump, any other type of electro-mechanical device, or brain shunts that contain magnets, keep sufficient distance (about 6 inches) between the device and the lamp.

Your lenses will yellow over time and should be replaced annually. Contact Verilux customer service at 1-800-786-6850 or order online at www.verilux.com.

Care and Cleaning

Your lamp is made from high-quality materials that will last for many years with minimum care. You may want to periodically clean the lamp using a mild non-abrasive cleaner and soft cloth. When cleaning, make sure you have turned off and unplugged the unit and allowed sufficient time for the unit to cool.

▲WARNING: To prevent risk of shock or personal injury when cleaning lamp, make sure you have turned off and unplugged the lamp and allowed sufficient time for the bulb to cool.

▲CAUTION: Clean only with a soft dry cloth – do not use solvents or cleaners containing abrasives or ammonia.

▲CAUTION: After cleaning the lamp, wipe away all moisture and let the lamp dry completely before restoring power.

Troubleshooting

Before requesting service on your Verilux lamp, please:

- Make sure all power connectors, cords and plugs are inserted fully and securely.
- Be sure the bulb is securely seated in the socket. (See pages 8 and 9.)
- Make sure there is power to the wall outlet, or try another outlet – one that is not connected to a dimmer.

PROBLEM	CHECK	SOLUTION
Light Won't Come On	Has the plug come loose?	Insert plug securely into wall outlet.
	Is the bulb old?	Check bulb for blackening/darkening at bottom of spiral where it joins the base (an indicator of bulb life) and replace bulb.
	Is the bulb connected securely into socket?	Be sure the bulb is securely seated into socket. Repeat the bulb installation process. DISCONNECT LAMP FROM OUTLET PRIOR TO REPLACING THE BULB IN INSTALLATION PROCESS.
Blackening/Darkening at Ends of Bulb	Is the bulb new?	Sometimes when a new bulb is turned on for the first time the tube will darken, but this will fade after a few minutes.
	Is the bulb old?	With older bulbs, darkening of the tube indicates end of bulb life. (This
Light Flickers	Has the plug come loose?	
	Is the bulb connected securely to the socket?	
	Is the lamp connected to a dimmer switch?	

HappyLight Liberty Series Natural Spectrum Energy Lamp

Environment: Operating Temperature: 10 to 30°C (50 to 86°F) Storage Temperature: -20 to 50°C (-4 to 122°F) Relative Humidity: 20 to 90% non-condensing

Bulb: Replacement bulb: 36 watt 10,000-hour bulb. Replace with Verilux brand only, Model number CFML36VLX by calling 1-800-454-4408 or visiting www.verilux.com.

Power: Input Power: 120 VAC, 60Hz, 0.50 A

ETL Listed in Canada and the USA

One-Year Limited Warranty

ATTENTION! ONCE OPENED, PLEASE DO NOT RETURN THIS PRODUCT TO THE STORE WHERE IT WAS PURCHASED FOR REPAIR OR REPLACEMENT!

Many questions can be answered by visiting www.verilux.com, or you may call our Customer Service Department at 1-800-786-6850 during normal business hours.

This limited warranty is provided by Verilux, Inc., 340 Mad River Park, Waitsfield, VT 05673

Verilux warrants this product to be free from defects in material and workmanship for a period of one year from the date of the original retail purchase from Verilux or an authorized Verilux distributor. **Proof of purchase is required for all warranty claims.** During the limited warranty period, Verilux, Inc., will, at its option, repair or replace defective parts of this product at no charge to the customer, subject to these limitations: This limited warranty does not include any postage, freight, handling, insurance or delivery fees. This warranty does not cover damage, defect or failure caused by or resulting from accident, external destruction, alteration, modification, abuse, misuse or misapplication of this product.

This warranty does not cover damage to the product resulting from return shipping or handling. Shipping insurance is recommended to help protect your product.

Return Authorization is required for all returns. To obtain a Return Authorization, please contact the Verilux Customer Service Department at 1-800-786-6850.

If, during the first year of ownership, this product fails to operate properly, it should be returned as specified at www.verilux.com/warrantyreplacement or as instructed by a Verilux customer service representative at 1-800-786-6850.

Note: *Verilux recommends using a quality surge suppressor on all electronic equipment. Voltage variations and spikes can damage electronic components in any system. A quality suppressor can eliminate the vast majority of failures attributed to surges. Suppressors may be purchased at electronics stores.*

Due to ongoing improvements, actual product may have slight variations from the product described in this manual.

Please visit our website at: www.verilux.com or call our Customer Service Department at
1-800-786-6850

Representatives are available Monday — Friday, 9am — 5pm EST

verilux®

340 Mad River Park, Waitsfield, VT 05673

Printed in China for: Verilux, Inc.

© Copyright 2017 Verilux, Inc. All rights reserved.

HappyLight®

Modèle : VT20

Edit#: 002	Title: VT20 Manual Trim to A5
Date: 09/06/17	Version: Rev20

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Esta garantía limitada es proporcionada por: Verilux, Inc., 340 Mad River Park, Waitsfield, VT 05673
Impreso en China para Verilux, Inc.